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August 21, 2026

Contents

1	Basic operation	3
1.1	Screens and settings	3
1.1.1	Main screen	3
1.1.2	Jog screen	5
1.1.3	Profile screen	7
1.1.4	Optical alignment screen	8
2	Clases de la API	12
2.1	Clase System	12
2.2	Clase AxisComponents	13
2.3	Clase Axis2D	14
2.4	Clase Axis3D	15
2.5	Clase Alignment	17
2.6	Clase Profile	20
2.7	Clase Angle Asjustment	21
3	Routines	23
3.1	RasterScan	23
3.1.1	clean_profile_data()	24
3.1.2	save_profile_txt()	25
3.1.3	create_heatmap()	26
3.2	Automapper	27
3.2.1	checkPositionsAutomap	27
3.2.2	moveAutomapOPM	28
3.2.3	Workflow	28
3.2.4	Relevant Configuration	29
3.2.5	Results	29
3.2.6	Precautions	29
3.2.7	moveAutomapOPM_OSA	29

3.2.8	Workflow	29
3.2.9	Relevant Configuration	30
3.2.10	Results	30
3.2.11	Precautions	30
A	Automap	31
B	Instalación TwinCAT XAE eXtended Automation Engineering	33
C	Documentacion de la API Suruga Seiki	39
C.1	Clase Alignment	39
C.2	Clase AngleAdjustment	46
C.3	Clase Axis2D	47
C.4	Clase Axis3D	48
C.5	Clase AxisComponents	49
C.6	Clase ErrorLog	49
C.7	Clase IO	50
C.8	Clase PowerMeter	50
C.9	Clase System	51

Chapter 1

Basic operation

1.1 Screens and settings

1.1.1 Main screen

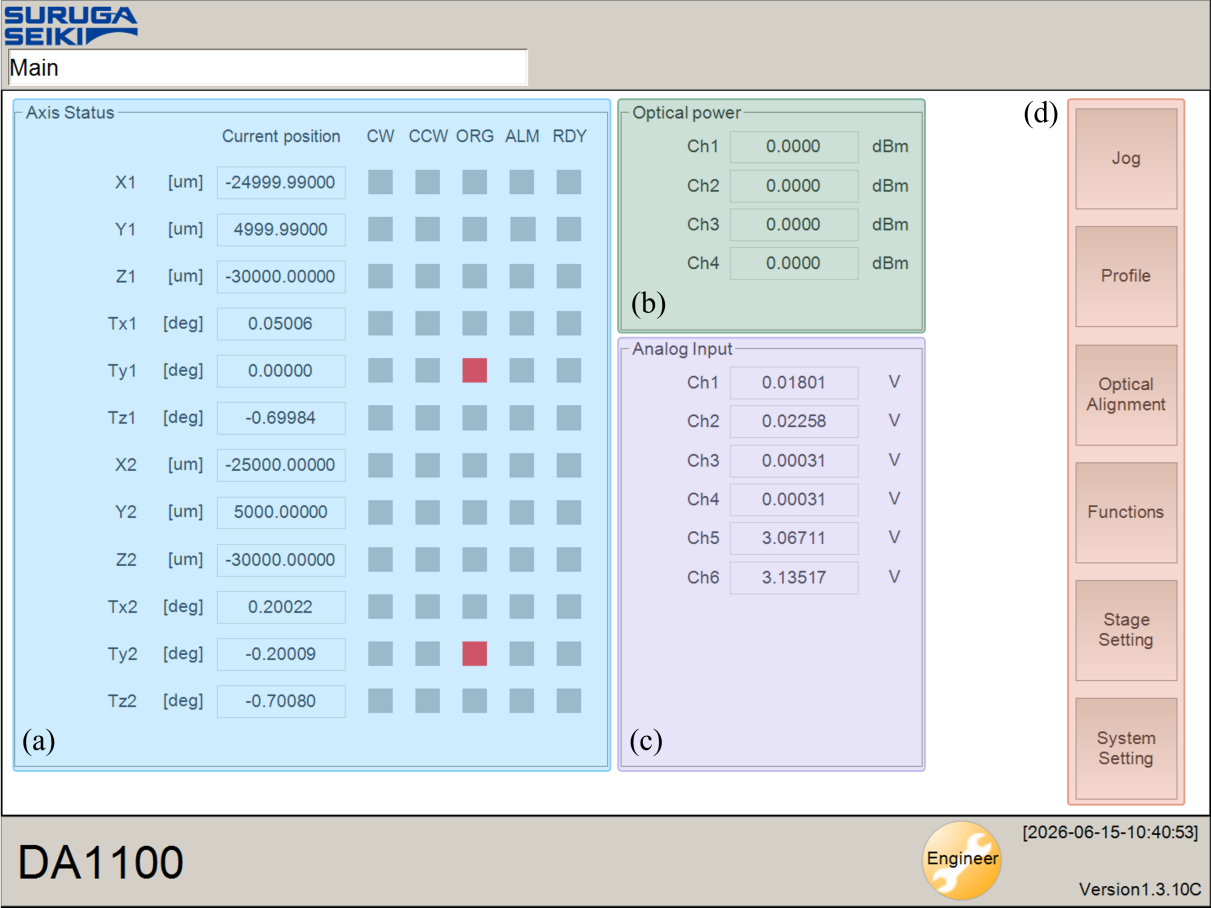
Figure 1.1: Main screen of the Suruga positioning system. The interface is divided into four main sections: (a) Axis Status, which displays the current position and status of each axis; (b) Optical Power, which displays the optical power measured by the connected power meter; (c) Analog Input, which displays the voltage measured at each available analog input channel; and (d) Function Access, which provides direct access to the main functionalities and configuration options of the Suruga platform.

This screen is divided into three main sections. The first section, Axis Status, allows the user to monitor the current position of each axis, the CW and CCW limit switch status, the origin sensor status, the alarm outputs, and the availability of each servo. The second section, Optical Power, displays the current reading of the power meter, if one is available. Finally, the Analog Input section allows the user to monitor the values of the different analog inputs. On the left side of the screen, direct access to the different functionalities of the Suruga platform is provided.

- **Current Position:** Displays the current coordinates of each axis. Linear axes are displayed in micrometers (μm), while rotational axes are displayed in degrees ($^\circ$). Axes with the number 1 correspond to the left arm of the Suruga positioning system, while axes with the number 2 correspond to the right arm.
- **CW:** Displays the status of the clockwise (CW) limit sensor. Red indicates that the limit has been reached, while gray indicates that the axis is within the safe operating range.
- **CCW:** Displays the status of the counterclockwise (CCW) limit sensor. Red indicates that the limit has been reached, while gray indicates that the axis is within the safe operating range.
- **ORG:** Displays the status of the origin sensor. For this particular stage, no origin sensor is available, so the indicator will remain continuously displayed in red or gray.
- **ALM:** Indicates the alarm status of the axis. Red is displayed when an alarm occurs, such as an overload or a similar fault condition. Gray indicates normal operation.
- **RDY:** Displays the operation status of each axis. Green indicates that the servo is enabled and ready for operation, while gray indicates that the servo is turned off.

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Routines

3.1 RasterScan

The routine performs a two-dimensional raster scan by moving along a selected axis while acquiring a profile along the perpendicular axis at each scan position. For example, when the X axis is selected as the scan axis, the system moves along X according to the specified scan step and acquires a Z-profile at each X position. This approach is particularly useful for vertical coupling measurements, where the optical power distribution can be characterized as a function of the X and Z positions.

During the scan, the profiles acquired using the Suruga Seiki utilities can be saved for further analysis. The acquired data are then cleaned by removing invalid or irrelevant samples that could interfere with the generation of the heatmap. The processed profiles are interpolated onto a common grid and used to generate an interactive heatmap representing the measured power distribution.

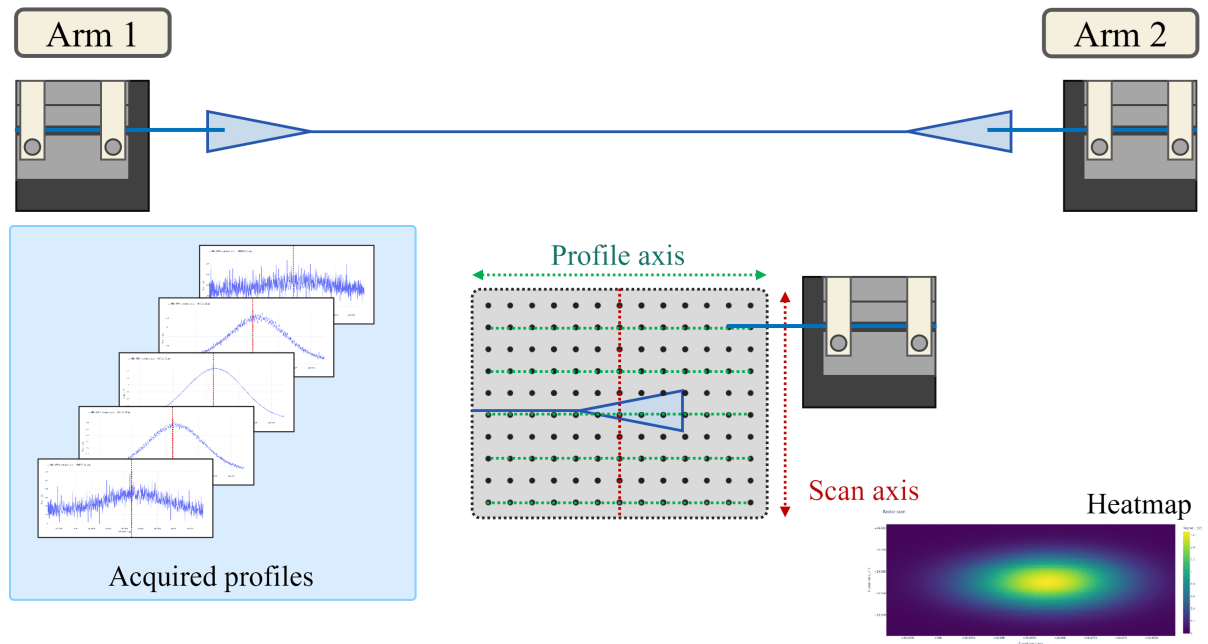
Figure 3.1: Two-dimensional raster scan for optical coupling localization. Profiles are acquired along the profile axis while the stage is incrementally moved along the perpendicular scan axis. The acquired profiles are combined to reconstruct a two-dimensional heatmap, allowing the position of maximum optical coupling to be identified.

Once the scan is complete, the routine identifies the global maximum power measured within the scanned area and, if enabled, automatically moves the arm to the corresponding X and Z position.

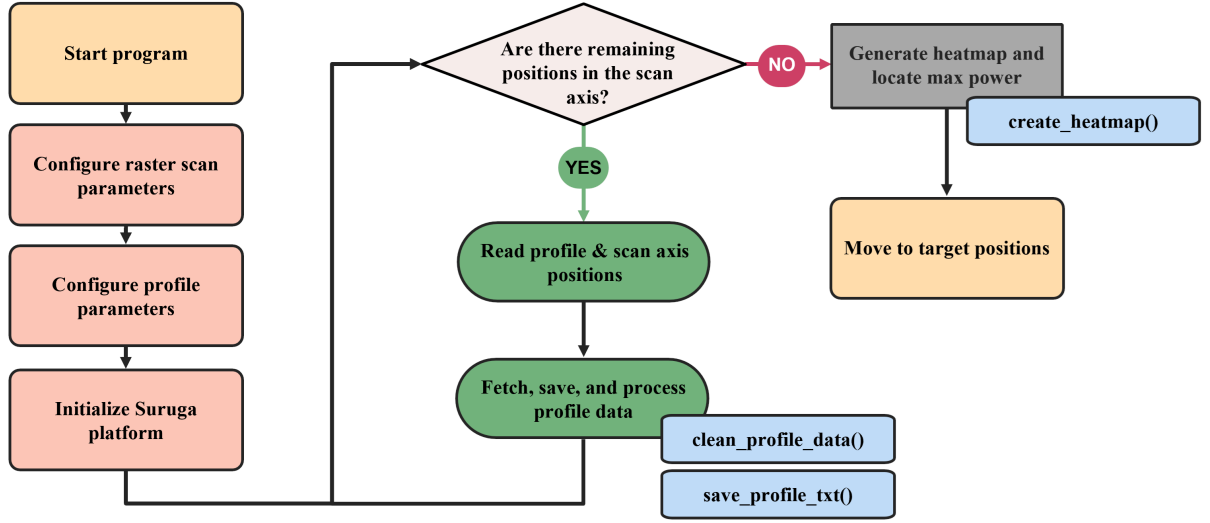
Figure 3.2: Flowchart of the *RasterScan* routine. The process begins with the configuration of the raster scan and profile parameters, followed by the initialization of the Suruga Seiki platform. For each position along the scan axis, the profile and scan-axis positions are obtained, and the acquired data are retrieved, stored, and processed using `clean_profile_data()` and `save_profile_txt()`. The process is repeated until all positions along the scan axis have been completed. Finally, the processed profiles are passed to `create_heatmap()` to reconstruct the two-dimensional power map and locate the maximum signal, after which the system can move to the corresponding target positions.

3.1.1 `clean_profile_data()`

Objective:

Clean and prepare the data acquired by the profile function for further processing.

This function receives the complete dictionary returned by `profile.get_profile_data()` and extracts only the variables required for the raster scan:

- **main_position**: the position of the axis being profiled.
- **signal1**: the measured optical signal.

The function then performs the following operations:

1. Converts the extracted data into NumPy arrays.
2. Removes invalid samples, including:
 - Non-finite position values.
 - Non-finite signal values.
 - Samples where the position is equal to zero.
3. Sorts and reorganizes the remaining samples according to their position.

The function **does not modify or save the original profile data**. Its purpose is to generate a cleaner and more manageable representation of the acquired data, which can then be used for further processing, particularly for generating the heatmap with `create_heatmap()`.

3.1.2 save_profile_txt()

Objective:

Preserve the complete profile data acquired during the measurement, together with the experimental parameters required to document and reproduce the measurement.

This function receives the complete dictionary returned by `profile.get_profile_data()`, including the following variables:

- **main_position:** position of the main axis during the profile.
- **sub1_position:** position of the first secondary axis.
- **sub2_position:** position of the second secondary axis.
- **signal1:** first measured signal.
- **signal2:** second measured signal.

In addition to the acquired data, the function constructs a header containing the experimental parameters and the stage positions at the time of acquisition:

- **Date:** Date and time of the measurement.
- **Center X:** X-axis center position.
- **Center Z:** Z-axis center position.
- **Arm:** Selected measurement arm or configuration.
- **Scan axis:** Axis along which the raster scan is performed.
- **Scan width:** Total width of the scan.
- **Scan step:** Distance between consecutive scan positions.
- **Profile axis:** Axis along which the profile is acquired.
- **Profile range:** Total range covered by the profile.
- **Profile speed:** Speed used during profile acquisition.

The function then stores the complete dataset and the associated metadata in a **tab-separated TXT file** using NumPy's `np.savetxt()` function.

Unlike `clean_profile_data()`, this function does **not** remove, filter, or reorganize the acquired samples. Its purpose is to preserve the original measurement data as completely as possible, while providing the experimental metadata necessary to identify and reproduce the measurement.

The resulting TXT file therefore contains both the **raw profile data** and the **measurement parameters** associated with its acquisition.

3.1.3 create_heatmap()

Objective:

Reconstruct a two-dimensional signal map from all the profiles acquired during the raster scan, generate an interactive heatmap, and determine the global maximum of `signal1` together with its corresponding X and Z coordinates.

The function receives two inputs:

- `scan_data`: a dictionary containing the processed profiles acquired during the raster scan. Each entry has the following structure:

```
[
  {
    "x": ...,
    "z_positions": ...,
    "signal": ...
  },
  ...
]
```

- `output_path`: the path where the generated heatmap will be saved.

Data Reconstruction

First, the function takes the Z positions from the first acquired profile and uses them as the common reference axis:

```
z_common = scan_data[0]["z_positions"]
```

Since the Z positions may not be identical for every profile, each profile is interpolated onto this common Z-axis. This ensures that all profiles have the same number of samples and can be represented consistently in a two-dimensional matrix.

For each profile, the function:

1. Retrieves the corresponding X position.
2. Retrieves the Z positions and measured signal.
3. Interpolates the signal onto `z_common`.
4. Adds the interpolated signal to the 2D signal matrix.

The resulting matrix represents the measured signal as a function of the X and Z coordinates and is used to generate the heatmap.

Heatmap Generation

The reconstructed 2D signal matrix is used to generate an **interactive Plotly heatmap**, where:

- The **X-axis** represents the raster scan position.
- The **Z-axis** represents the profile position.
- The **color scale** represents the measured `signal1`.

The resulting heatmap provides a two-dimensional visualization of the measured optical signal across the scanned area and is saved to the location specified by `output_path`.

Global Maximum Detection

In addition to generating the heatmap, the function determines the position of the maximum measured signal.

For each profile, the function first identifies its **local maximum**. The local maximum is then compared with the maximum found in the previous profiles. By continuously retaining the highest value, the function ultimately determines the **global maximum** across the entire raster scan.

The global maximum is stored together with its corresponding coordinates:

- **Maximum signal:** maximum value of `signal1`.
- **X coordinate:** X position of the profile where the maximum was found.
- **Z coordinate:** Z position corresponding to the maximum signal within that profile.

Therefore, `create_heatmap()` performs two main tasks: it **reconstructs and visualizes the complete raster scan as a two-dimensional signal map**, and it **identifies the position of the maximum output power** within the scanned region.

3.2 Automapper

3.2.1 checkPositionsAutomap

Validates the measurement plan generated by AutoMapper before moving the Suruga platform. It calculates and displays the absolute `x1` and `x2` positions for each device defined in the Excel file, but it does not connect to any instrument or move any axes.

Source code: [checkPositionsAutomap.py](#).

Workflow

1. Loads the chip layout and measurement plan from an Excel file using `automaper.load_chip()`.
2. Uses `origin_x1` and `origin_x2` as the reference origins.
3. For each measurement, calculates the input and output offsets and the absolute positions of both arms.
4. Displays the device, ports, and calculated positions, and waits for confirmation before displaying the next measurement.

Relevant Configuration

Table 3.1: Relevant configuration parameters for `checkPositionsAutomap`.

Parameter	Initial Value	Description
EXCEL_FILE	chipMaps/SurugaTest.xlsx	Chip layout and AutoMapper measurement plan.
origin_x1 / origin_x2	0.0 μm / 0.0 μm	Reference positions used to calculate the absolute coordinates.

The absolute coordinates are calculated according to

$$x_1 = \text{origin}_{x1} - \text{input}.x \quad (3.1)$$

and

$$x_2 = \text{origin}_{x2} - \text{output}.x. \quad (3.2)$$

Recommended Usage

Run this routine whenever the Excel file, chip coupling, or origin convention changes. Compare the printed positions against the mechanical setup limits before running a routine that moves the machine.

Results

This routine does not create files or perform measurements. Its output consists of the position list displayed in the console.

3.2.2 moveAutomapOPM

Automatically iterates through the measurements defined by AutoMapper, positions the fibers at the input and output waveguides, performs alignment on both arms, and measures the average optical power using the optical power meter (OPM).

Source code: [moveAutomapOPM.py](#).

3.2.3 Workflow

1. Loads the measurement plan from the AutoMapper Excel file.
2. Connects to the Suruga platform and creates controllers for `x1`, `x2`, alignment, and the power meter.
3. Loads `origin.txt` or uses the current position as the WGO origin.
4. For each device, calculates the absolute positions of both arms and moves `x1` and `x2`.
5. Performs a *flat* alignment on the left arm followed by another alignment on the right arm.
6. Reads ten samples from the OPM, calculates their average, and displays the result in the console.

7. Requests confirmation before continuing and, once all measurements are completed, returns to WG0.

3.2.4 Relevant Configuration

Table 3.2: Relevant configuration parameters for `moveAutomapOPM`.

Parameter	Initial Value	Description
EXCEL_FILE	chipMaps/SurugaTest.xlsx	Chip layout and measurement plan.
USE_SAVED_ORIGIN / SAVE_ORIGIN	True / False	Control whether <code>origin.txt</code> is loaded or saved.
ORIGIN_FILE	origin.txt	File containing the X1 and X2 coordinates of WG0.
flat_left / flat_right	---	Alignment parameters for each arm.

3.2.5 Results

The average optical power for each device is printed to the console. The routine does not save the measurement values to a file.

3.2.6 Precautions

The routine moves both arms and requires manual confirmation between measurements. Verify the position plan using `checkPositionsAutomap` before running it.

3.2.7 `moveAutomapOPM_OSA`

Characterizes the devices defined by AutoMapper by combining Suruga platform movement, automatic alignment, an average optical power measurement using the OPM, and spectrum acquisition with the OSA.

Source code: `moveAutomapOPM_OSA.py`.

3.2.8 Workflow

1. Loads the chip layout and measurement plan from Excel.
2. Connects to the Suruga platform, OPM, and OSA, and configures the spectral sweep.
3. Loads or saves the WG0 origin in `origin2.txt`.
4. For each device and repetition, positions x1 and x2 based on the corresponding input and output waveguides.
5. Aligns the left and right arms.
6. Runs an OSA sweep, retrieves the trace, and averages ten power readings.
7. Saves the trace along with the device, ports, optical power, and repetition number as metadata.
8. Returns to WG0 and closes the OSA connection.

3.2.9 Relevant Configuration

Table 3.3: Relevant configuration parameters for `moveAutomapOPM_OSA`.

Parameter	Initial Value	Description
DEBUG_MODE	True	Requests manual confirmations; when set to <code>False</code> , the routine runs automatically.
USE_SAVED_ORIGIN / SAVE_ORIGIN	False / True	Control whether the origin stored in <code>origin2.txt</code> is loaded or saved.
EXCEL_FILE	chipMaps/SurugaTest01.xlsx	Chip layout and AutoMapper measurement plan.
num.repeticiones	1	Number of complete passes through the measurement plan.
CENTER / SPAN	1550 nm / 90 nm	OSA spectral window.
RESOLUTION / SENSITIVITY	0.2 nm / -65 dBm	OSA resolution and sensitivity.
channel	1	Power meter channel.

3.2.10 Results

For each device, the routine saves a trace to:

```
data/<date>/medicion_<repetition>_<device>_<input>_<output>.txt
```

The file contains the spectral data and metadata including the device, ports, average optical power, and repetition number.

3.2.11 Precautions

With `DEBUG_MODE=True`, the routine pauses before the OSA sweep, between devices, and before returning to the origin.

Before using automatic mode, validate the movements with `checkPositionsAutomap`, verify the origin, and confirm that the OSA VISA resource is available.

Appendix A

Automap

Automaper is a python script that converts a simple Excel description of a photonic chip into a complete measurement plan ready for chip characterization.

Instead of hardcoding device coordinates or manually programming measurement routines, the user only needs to describe the chip structure in an Excel file then automaper.py interprets this information and automatically generates the data structures required by the positioning software.

The Excel file is organized in independent sheets, where each sheet describes a different aspect of the chip as follows:

Metadata

Contains general information about the chip, answers the question: What is the context of this measurement?.

Typical information includes

- Project
- Wafer
- Die
- Operator
- Date
- Comments

This information is descriptive only and is not involved in the positioning process.

Devices

This worksheet answers the question: What devices are available on the chip for measurement?

Each row defines a single device together with the information required to identify it throughout the measurement plan.

The columns have the following purposes:

- **Device:** A unique identifier for each device. This identifier is used throughout the workbook to associate the corresponding ports and measurement paths with the correct device.

- **Type:** Specifies the type of photonic device (e.g., waveguide, MMI, ring resonator, splitter). Although this information is not required for positioning, it is useful for filtering devices, generating reports, and applying device-specific measurement or analysis routines.
- **Description:** An optional user-defined description. This field does not affect Automaper's operation or the alignment process; it is provided solely for documentation and traceability purposes.
- **Measure:** Enables or disables a device for the current measurement routine. When unchecked, the device is excluded from the measurement routine without requiring any modification to the measurement paths.

Ports

Where is every optical port located?

Measurement plan

Which optical paths should be measured?

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C.2 Clase AngleAdjustment

```

=====
SurugaSeiki.Motion.AngleAdjustment
=====

```

```

Constructores:
    AngleAdjustment(UInt16 id)

Metodos:
    ErrorCode Start()
    ErrorCode Stop()
    Status GetStatus()
    UInt32 GetErrorAxisID()
    AdjustingStatus GetAdjustingStatus()
    UInt32 GetProfilePacketSumIndex(ProfileDataType profileDataType)
    ProfileData RequestProfileData(ProfileDataType profileDataType, UInt32 index)
    ErrorCode SetParameter(AngleAdjustmentParameter parameter)
    ErrorCode SetParameter(AngleAdjustmentParameter parameter, RotationCenter rotationCenter)

```

```

Tipos anidados/enums:
    SurugaSeiki.Motion.AngleAdjustment+ErrorCode
    SurugaSeiki.Motion.AngleAdjustment+Status
    SurugaSeiki.Motion.AngleAdjustment+AdjustingStatus
    SurugaSeiki.Motion.AngleAdjustment+ProfileDataType
    SurugaSeiki.Motion.AngleAdjustment+AngleAdjustmentParameter
    SurugaSeiki.Motion.AngleAdjustment+ProfileData

```

```

=====
SurugaSeiki.Motion.AngleAdjustment+AdjustingStatus
=====

```

```

Campos:
    Int32 value_
    static AdjustingStatus NotAdjusting
    static AdjustingStatus Initializing
    static AdjustingStatus ContactingZ
    static AdjustingStatus AdjustingTx
    static AdjustingStatus AdjustingTy

```

```

=====
SurugaSeiki.Motion.AngleAdjustment+AngleAdjustmentParameter
=====

```

```

Campos:
    Double gap
    UInt16 signalChNumber

```

```

Double signalLowerLimit
UInt16 unlockDOutChNumber
Boolean unlockDOutControlOn
Boolean lockUnlockAdjustEnbale
Double lockUnlockDifference
UInt16 contactAxisNumber
Double contactSearchRange
Double contactSearchSpeed
UInt32 contactSmoothing
UInt32 contactSensitivity
Double pushDistance
UInt16 angleAxisNumberTx
UInt16 angleAxisNumberTy
Double angleSearchRangeTx
Double angleSearchRangeTy
Double angleSearchSpeedTx
Double angleSearchSpeedTy
UInt32 angleSmoothingTx
UInt32 angleSmoothingTy
UInt32 angleSensitivityTx
UInt32 angleSensitivityTy
UInt16 angleJudgeCountTx
UInt16 angleJudgeCountTy
Double angleConvergentRangeTx
Double angleConvergentRangeTy
UInt16 angleComparisonCount
UInt16 angleMaxCount
=====
SurugaSeiki.Motion.AngleAdjustment+ErrorCode
=====

Campos:
    Int32 value__
    static ErrorCode None
    static ErrorCode Axis
    static ErrorCode Adjusting
    static ErrorCode Parameter
=====
SurugaSeiki.Motion.AngleAdjustment+ProfileData
=====

Campos:
    UInt16 mainAxisNumber
    UInt16 sub1AxisNumber
    UInt16 sub2AxisNumber
    UInt16 signalCh1Number
    UInt16 signalCh2Number
    UInt32 packetIndex
    UInt32 dataCount
    Double[] mainPositionList
    Double[] sub1PositionList
    Double[] sub2PositionList
    Double[] signalCh1List
    Double[] signalCh2List
=====
SurugaSeiki.Motion.AngleAdjustment+ProfileDataType
=====

Campos:
    Int32 value__
    static ProfileDataType ContactZ
    static ProfileDataType AdjustmentTx
    static ProfileDataType AdjustmentTy
=====
SurugaSeiki.Motion.AngleAdjustment+Status
=====

Campos:
    Int32 value__
    static Status Stopping
    static Status Success
    static Status Adjusting
    static Status ProfileDataOver
    static Status InvalidParameter
    static Status ServoIsNotReady
    static Status ServoIsAlarm
    static Status StageOnLimit
    static Status SignalLowerLimit
    static Status CouldnotContact
    static Status AdjustCountOver
    static Status AngleAdjustRnageOver
    static Status LostContact

```

C.3 Clase Axis2D

```

=====
SurugaSeiki.Motion.Axis2D
=====

```

```
Constructores:
    Axis2D(UInt16 mainStageNumber, UInt16 subStageNumber)

Metodos:
    ErrorCode SetAxisNumber(UInt16 mainStageNumber, UInt16 subStageNumber)
    UInt16[] GetAxisNumber()
    ErrorCode MoveAbsolute(Point targetPosition)
    ErrorCode MoveAbsolute(Double mainTargetPosition, Double subTargetPosition)
    ErrorCode MoveRelative(Point targetDistance)
    ErrorCode MoveRelative(Double mainDistance, Double subDistance)
    ErrorCode MoveAbsoluteAngle(Double mainTargetPosition, Double subAngle)
    ErrorCode MoveRelativeAngle(Double mainDistance, Double subAngle)
    ErrorCode Stop()
    Point GetTargetPosition()
    Point GetActualPosition()
    Boolean IsMoving()
    ErrorCode SetSpeed(Double speed)
    ErrorCode SetAccelRate(Double rate)
    ErrorCode SetDecelRate(Double rate)
    ErrorCode GetErrorCode()
    Status GetStatus()
```

```
Tipos anidados/enums:
    SurugaSeiki.Motion.Axis2D+Point
```

```
=====
SurugaSeiki.Motion.Axis2D+Point
=====
```

```
Campos:
    Double X
    Double Y
```

C.4 Clase Axis3D

```
=====
SurugaSeiki.Motion.Axis3D
=====
```

```
Constructores:
    Axis3D(UInt16 mainStageNumberX, UInt16 mainStageNumberY, UInt16 mainStageNumberZ)
```

```
Metodos:
    ErrorCode SetAxisNumber(UInt16 mainStageNumberX, UInt16 mainStageNumberY, UInt16 mainStageNumberZ)
    UInt16[] GetAxisNumber()
    ErrorCode MoveAbsolute(Point3D targetPosition)
    ErrorCode MoveAbsolute(Double targetPositionX, Double targetPositionY, Double targetPositionZ)
    ErrorCode MoveAbsolute(Point3D targetPosition, UInt16 dOutId, Int16 dOutOnDelayTime, Int16 dOutOffDelayTime)
    ErrorCode MoveAbsolute(Double targetPositionX, Double targetPositionY, Double targetPositionZ, UInt16 dOutId, Int16 dOutOnDelayTime, Int16 dOutOffDelayTime)
    ErrorCode MoveRelative(Point3D targetDistance)
    ErrorCode MoveRelative(Double targetDistanceX, Double targetDistanceY, Double targetDistanceZ)
    ErrorCode MoveRelative(Point3D targetDistance, UInt16 dOutId, Int16 dOutOnDelayTime, Int16 dOutOffDelayTime)
    ErrorCode MoveRelative(Double targetDistanceX, Double targetDistanceY, Double targetDistanceZ, UInt16 dOutId, Int16 dOutOnDelayTime, Int16 dOutOffDelayTime)
    ErrorCode Stop()
    Point3D GetTargetPosition()
    Point3D GetActualPosition()
    Boolean IsMoving()
    ErrorCode SetSpeed(Double speed)
    ErrorCode SetAccelRate(Double accelRate)
    ErrorCode SetDecelRate(Double decelRate)
    ErrorCode SetStageNumber(UInt16 mainStageNumberX, UInt16 mainStageNumberY, UInt16 mainStageNumberZ)
    ErrorCode SetSubStageNumber(UInt16 subStageNumberX, UInt16 subStageNumberY, UInt16 subStageNumberZ)
    ErrorCode SetRotationCenter(Double rotationOffsetX, Double rotationOffsetY, Double rotationOffsetZ)
    ErrorCode GetErrorCode()
    Status GetStatus()
```

```
Tipos anidados/enums:
    SurugaSeiki.Motion.Axis3D+Point3D
    SurugaSeiki.Motion.Axis3D+RotationCenter
```

```
=====
SurugaSeiki.Motion.Axis3D+Point3D
=====
```

```
Campos:
    Double X
    Double Y
    Double Z
```

```
=====
SurugaSeiki.Motion.Axis3D+RotationCenter
=====
```

```
Campos:
    UInt16 mainStageNumberX
    UInt16 mainStageNumberY
    UInt16 mainStageNumberZ
    UInt16 subStageNumberX
    UInt16 subStageNumberY
    UInt16 subStageNumberZ
    Double offsetX
```

```
Double offsetY
Double offsetZ
```

C.5 Clase AxisComponents

```
=====
SurugaSeiki.Motion.AxisComponents
=====

Constructores:
    AxisComponents(UInt16 axisNumbern)

Metodos:
    ErrorCode MoveAbsolute(Double position)
    ErrorCode MoveRelative(Double distance)
    ErrorCode ReturnOrigin()
    ErrorCode Stop()
    ErrorCode TurnOnServo()
    ErrorCode TurnOffServo()
    ErrorCode ResetError()
    Double GetTargetPosition()
    Double GetActualPosition()
    Boolean IsServoOn()
    Boolean IsServoAlarm()
    Boolean IsMoving()
    Boolean IsOriginSignalOn()
    Boolean IsCwLimitSignalOn()
    Boolean IsCcwLimitSignalOn()
    Boolean IsCwTorqueLimitOn()
    Boolean IsCcwTorqueLimitOn()
    Boolean IsOriginReturnSucceeded()
    ErrorCode SetAccelRate(Double rate)
    ErrorCode SetDecelRate(Double rate)
    ErrorCode SetSineMotion(Boolean sineMotion)
    ErrorCode SetFulcrumSpan(Double span)
    ErrorCode SetEncoderStepLength(Double length)
    ErrorCode SetMaxSpeed(Double speed)
    ErrorCode SetCWTorqueLimit(Double cWlimit)
    ErrorCode SetCCWTorqueLimit(Double cCWlimit)
    ErrorCode GetErrorCode()
    Status GetStatus()

Tipos anidados/enums:
    SurugaSeiki.Motion.AxisComponents+ErrorCode
    SurugaSeiki.Motion.AxisComponents+Status
```

```
=====
SurugaSeiki.Motion.AxisComponents+ErrorCode
=====
```

```
Campos:
    Int32 value__
    static ErrorCode None
    static ErrorCode ServoAlarm
    static ErrorCode ServoOff
    static ErrorCode Moving
    static ErrorCode NotConnected
    static ErrorCode Limit
    static ErrorCode Origin
    static ErrorCode Parameter
    static ErrorCode Torque
```

```
=====
SurugaSeiki.Motion.AxisComponents+Status
=====
```

```
Campos:
    Int32 value__
    static Status ServoAlarm
    static Status ServoOff
    static Status Stopping
    static Status InPosition
    static Status Moving
    static Status NotConnected
    static Status CwLimit
    static Status CcwLimit
    static Status CwTorqueLimit
    static Status CcwTorqueLimit
```

C.6 Clase ErrorLog

```
=====
SurugaSeiki.Motion.ErrorLog
=====
```

C.7 Clase IO

```
=====
SurugaSeiki.Motion.IO
=====

Constructores:
    IO()

Metodos:
    Boolean GetPortState(DigitalIOType iOType, UInt16 id)
    Void SetPortState(UInt16 id, Boolean portState)
    Void SetPulseOutput(UInt16 id, UInt16 delayTime, UInt16 highWidth, UInt16 lowWidth, Int16 repetitionCount)
    Double GetAnalogValue(AnalogIOType iOType, UInt16 id)
    Void SetAnalogValue(UInt16 id, Double setValue)

Campos:
    static Int16 RepetitionCountForever

Tipos anidados/enums:
    SurugaSeiki.Motion.IO+DigitalIOType
    SurugaSeiki.Motion.IO+AnalogIOType

=====
SurugaSeiki.Motion.IO+AnalogIOType
=====

Campos:
    Int32 value_
    static AnalogIOType Input
    static AnalogIOType Output

=====
SurugaSeiki.Motion.IO+DOutControlParameterToAlignmentSystem
=====

=====
SurugaSeiki.Motion.IO+DigitalIOType
=====

Campos:
    Int32 value_
    static DigitalIOType Input
    static DigitalIOType Output
```

C.8 Clase PowerMeter

```
=====
SurugaSeiki.Motion.PowerMeter
=====

Constructores:
    PowerMeter()

Metodos:
    Double GetPower(UInt16 ch)
    AllowedRange GetRange(UInt16 ch)
    Double GetWavelength(UInt16 ch)
    ErrorCode SetRange(UInt16 ch, AllowedRange range)
    ErrorCode SetWaveLength(UInt16 ch, Double waveLength)
    ErrorCode PauseCommunication(UInt16 ch, Boolean isEnabled)

Tipos anidados/enums:
    SurugaSeiki.Motion.PowerMeter+ErrorCode
    SurugaSeiki.Motion.PowerMeter+AllowedRange

=====
SurugaSeiki.Motion.PowerMeter+AllowedRange
=====

Campos:
    Int16 value_
    static AllowedRange Minus60
    static AllowedRange Minus50
    static AllowedRange Minus40
    static AllowedRange Minus30
    static AllowedRange Minus20
    static AllowedRange Minus10
    static AllowedRange Zero
    static AllowedRange Plus10
    static AllowedRange Plus20
    static AllowedRange Auto
    static AllowedRange Invalid

=====
SurugaSeiki.Motion.PowerMeter+ErrorCode
=====

Campos:
    Int32 value_
    static ErrorCode None
    static ErrorCode ChannelOutOfRange
```

```
static ErrorCode Invalid
```

C.9 Clase System

```
=====
SurugaSeiki.Motion.System
=====
```

Propiedades:

```
UInt16 MaxAxisStageNumber (get/set)
UInt16 MaxAnalogInputChNumber (get/set)
UInt16 MaxAlignmentNumber (get/set)
UInt16 MaxAngleAdjustmentNumber (get/set)
UInt16 MaxDInChNumber (get/set)
UInt16 MaxDOutChNumber (get/set)
UInt16 MaxAnalogOutputChNumber (get/set)
UInt16 MaxPowerMeterChNumber (get/set)
System Instance (get)
Boolean Connected (get)
String DllVersion (get)
String SystemVersion (get)
```

Metodos:

```
Boolean IsError()
ErrorCode GetCurrentErrorCode()
Boolean IsEmergencyAsserted()
Void SetAddress(String addr)
Void DisconnectAds()
```

Campos:

```
static UInt16 MaxProfileNumber
static UInt32 ProfileDataPacketMaxDataNumber
```

Tipos anidados/enums:

```
SurugaSeiki.Motion.System+ErrorCode
```

```
=====
SurugaSeiki.Motion.System+AlignmentChParam
=====
```

```
=====
SurugaSeiki.Motion.System+AlignmentCommand
=====
```

Campos:

```
Int32 value__
static AlignmentCommand Null
static AlignmentCommand StartSingle
static AlignmentCommand StartFlat
static AlignmentCommand StartFocus
static AlignmentCommand StartTwoChRotation
static AlignmentCommand StartMMFRotation
static AlignmentCommand Stop
static AlignmentCommand SetSingle
static AlignmentCommand SetFlat
static AlignmentCommand SetFocus
static AlignmentCommand SetTwoChRotation
static AlignmentCommand SetMMFRotation
static AlignmentCommand StartLoopBack
static AlignmentCommand SetLoopBack
```

```
=====
SurugaSeiki.Motion.System+AlignmentMainAxisParam
=====
```

```
=====
SurugaSeiki.Motion.System+AlignmentParam
=====
```

```
=====
SurugaSeiki.Motion.System+AlignmentPowerMeterInfo
=====
```

```
=====
SurugaSeiki.Motion.System+AlignmentStatusInfo
=====
```

```
=====
SurugaSeiki.Motion.System+AnalogOutputControlCommand
=====
```

Campos:

```
Int32 value__
static AnalogOutputControlCommand Null
static AnalogOutputControlCommand SetAnalogOutput
```

```
=====
SurugaSeiki.Motion.System+AngleAdjustmentCommand
=====
```

Campos:

```
Int32 value__
static AngleAdjustmentCommand Null
```

```

static AngleAdjustmentCommand SetParameter
static AngleAdjustmentCommand Start
static AngleAdjustmentCommand Stop

=====
SurugaSeiki.Motion.System+AngleAdjustmentParam
=====

SurugaSeiki.Motion.System+AngleAdjustmentStatusInfo
=====

SurugaSeiki.Motion.System+Axis2DCommand
=====

Campos:
  Int32 value__
  static Axis2DCommand Null
  static Axis2DCommand MoveAbsolute
  static Axis2DCommand MoveRelative
  static Axis2DCommand MoveAbsoluteAngle
  static Axis2DCommand MoveRelativeAngle
  static Axis2DCommand Stop
  static Axis2DCommand SetCompositionSpeed
  static Axis2DCommand SetAccelRate
  static Axis2DCommand SetDecelRate
  static Axis2DCommand SetStageNumber

=====
SurugaSeiki.Motion.System+Axis2DParam
=====

SurugaSeiki.Motion.System+AxisCommandStruct
=====

SurugaSeiki.Motion.System+AxisComponentCommand
=====

Campos:
  Int32 value__
  static AxisComponentCommand Null
  static AxisComponentCommand MoveAbsolute
  static AxisComponentCommand MoveRelative
  static AxisComponentCommand ReturnOrigin
  static AxisComponentCommand Stop
  static AxisComponentCommand TurnOnServo
  static AxisComponentCommand TurnOffServo
  static AxisComponentCommand ResetError
  static AxisComponentCommand SetMaxSpeed
  static AxisComponentCommand SetAccelRate
  static AxisComponentCommand SetDecelRate
  static AxisComponentCommand SetSineMotion
  static AxisComponentCommand SetFulcrumSpan
  static AxisComponentCommand SetEncoderStepLength
  static AxisComponentCommand SetCWTorquelimit
  static AxisComponentCommand SetCCWTorquelimit

=====
SurugaSeiki.Motion.System+AxisControlType
=====

Campos:
  Int32 value__
  static AxisControlType AxisComponents
  static AxisControlType Axis2D
  static AxisControlType Axis3D

=====
SurugaSeiki.Motion.System+AxisMethod
=====

Campos:
  Int32 value__
  static AxisMethod Null
  static AxisMethod MoveAbsolute
  static AxisMethod MoveRelative
  static AxisMethod MoveAbsoluteAngle
  static AxisMethod MoveRelativeAngle
  static AxisMethod Stop
  static AxisMethod SetCompositionSpeed
  static AxisMethod SetAccelRate
  static AxisMethod SetDecelRate
  static AxisMethod SetStageNumber
  static AxisMethod SetSubStageNumber
  static AxisMethod SetRotationCenter
  static AxisMethod SetFixedAxisOffset

=====
SurugaSeiki.Motion.System+AxisStatus
=====

SurugaSeiki.Motion.System+DOutControlCommand
=====

Campos:

```

```

    Int32 value__
    static DOutControlCommand Null
    static DOutControlCommand SetDout
    static DOutControlCommand SetPulseOutput

=====
SurugaSeiki.Motion.System+DeviceControlCommand
=====

Campos:
    Int32 value__
    static DeviceControlCommand Null
    static DeviceControlCommand SetWaveLength

=====
SurugaSeiki.Motion.System+DeviceStatusMap
=====

=====
SurugaSeiki.Motion.System+ErrorCode
=====

Campos:
    Int32 value__
    static ErrorCode NoError
    static ErrorCode ConnectionLost
    static ErrorCode VersionError
    static ErrorCode MaxNumberOfErrorCode

=====
SurugaSeiki.Motion.System+Point2d
=====

=====
SurugaSeiki.Motion.System+PowerMeterControlCommand
=====

Campos:
    Int32 value__
    static PowerMeterControlCommand None
    static PowerMeterControlCommand SetRange
    static PowerMeterControlCommand SetWaveLength
    static PowerMeterControlCommand SetCommunicationPaused

=====
SurugaSeiki.Motion.System+ProfileCommand
=====

Campos:
    Int32 value__
    static ProfileCommand Null
    static ProfileCommand SetParameter
    static ProfileCommand Start
    static ProfileCommand Stop

=====
SurugaSeiki.Motion.System+ProfileParam
=====

=====
SurugaSeiki.Motion.System+RotationCenterParam
=====

=====
SurugaSeiki.Motion.System+SystemInfo
=====

```

Bibliography